

The logo features the text "SS SYSTEM" in a bold, orange, sans-serif font. The text is centered within a large, dark circular shape that has a glowing orange and yellow ring around its right edge, resembling a sunset or a bright light source. The background of the circle is black, and the overall effect is a high-contrast, modern design.

SS SYSTEM

IMPLANT

- 090** SSII SA Implant
- 092** SSIII SA Implant
- 096** Simple Mount
- 096** Cover Screw
- 097** Closing Screw
- 098** Healing Abutment


COMPONENTS

- 100** PROSTHETIC FLOW DIAGRAM 1
- 101** Solid Abutment
- 104** Excellent Solid Abutment
- 108** PROSTHETIC FLOW DIAGRAM 2
- 109** ComOcta Abutment
- 114** ComOcta Plus Abutment
- 116** ComOcta Milling Abutment
- 117** ComOcta Gold Abutment
- 118** ComOcta NP-Cast Abutment
- 119** ComOcta Temporary Abutment
- 120** ComOcta Angled Abutment
- 122** PROSTHETIC FLOW DIAGRAM 3
- 123** Octa Abutment
- 128** PROSTHETIC FLOW DIAGRAM 4
- 129** Port Abutment
- 132** Locator® Abutment
- 136** O-ring Abutment
- 138** OneSeal

SS Design & Surface Feature



SS

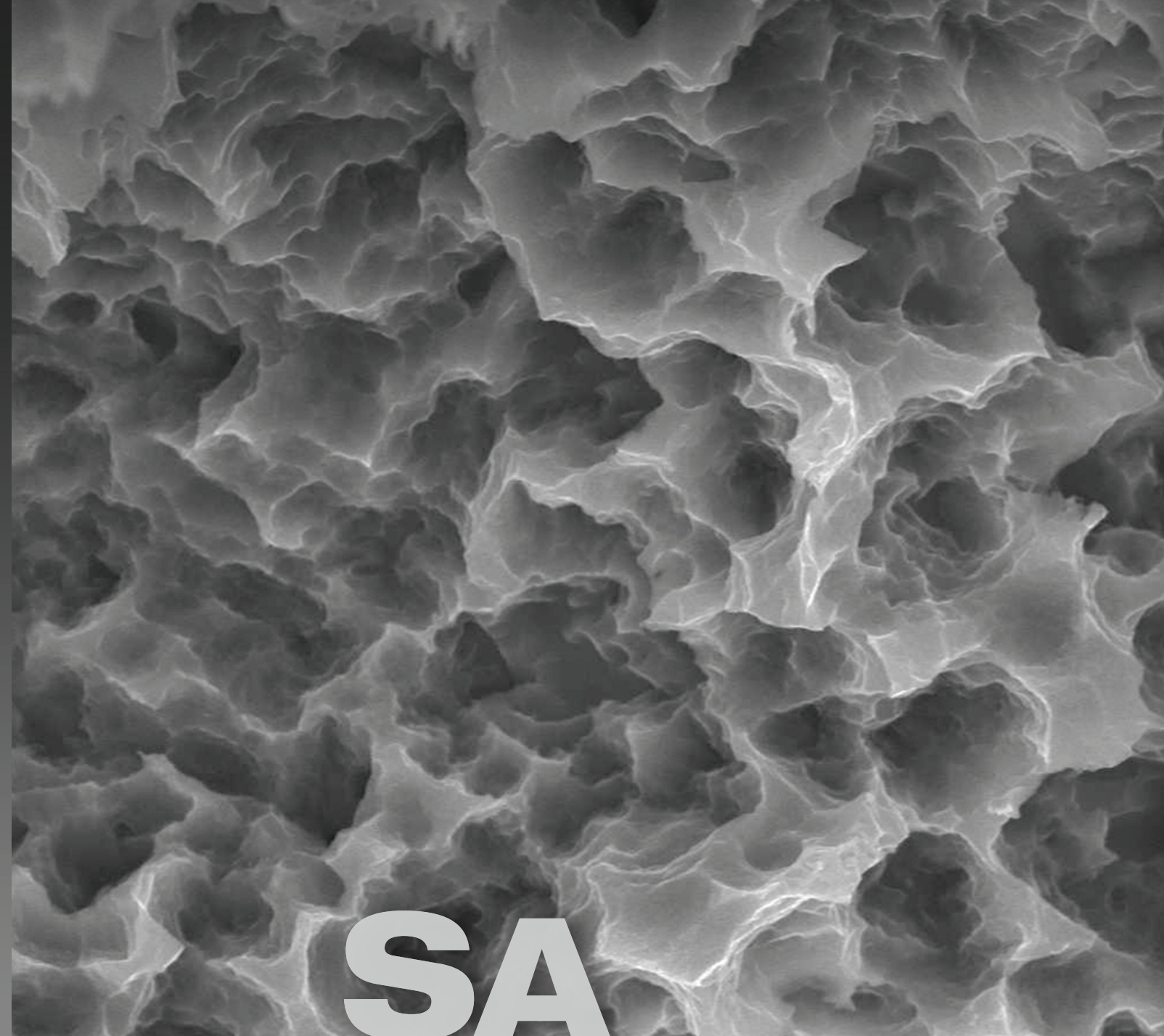

Non-submerged type implant with an internal octa 8° tapered connection based on the 1st stage surgery

- Connection - Regular / Wide
- Corkscrew thread & cutting edge
 - Superior self-threading effect for ease of placement path adjustment
 - Enhanced initial stability in soft bone and application of consistent placement torque according to the drill diameter
- Various body shape options available to match the patient's bone quality and clinical condition
 - SSII (straight body) : Ease of placement depth adjustment
 - SSIII (1.5° tapered body) : Excellent initial stability needed for immediate loading even in soft bone
- Applicable surface types - SA

Optimized surface morphology through acid-etching treatment

- Surface roughness: Ra 2.0-3.0 μ m (Note: the roughness in the upper 0.5mm part is Ra 0.5-0.6 μ m)
- Uniform surface micro-pits of 1~3 μ m
- Surface area increased by 46% compared to resorbable blast media (RBM) treated implants

In-vitro and In-vivo Bone Response

- Osteoblast differentiation and ossification improved by 20% compared to RBM-treated implants
- Initial bone response in a large animal model (mini-pig)
 - Initial stability (removal torque (RT), 4 weeks) improved by 48% compared to RBM-treated implants
 - Ossification (bone implant contact (BIC), 4 weeks) improved by 20% compared to RBM-treated implants



SS packaging color information

SSII SA Implant


- Non-submerged type implant with an internal octa 8° tapered connection based on the one-stage surgery
- Optimal thread design for realization of optimal SA surface
- Straight body design for easy adjustment of placement depth
- Superior self-threading effect with corkscrew thread
- Recommended placement torque : ≤ 40Ncm
- ※ Implants with a diameter of 4.5mm or greater are recommended for the posterior region with a single case

NoMount implant order code

: implant product code (ex : SS2R4011S18)

Pre-Mounted implant (implant + simple mount + cover screw) order code

: A + implant product code (ex : ASS2R4011S18)

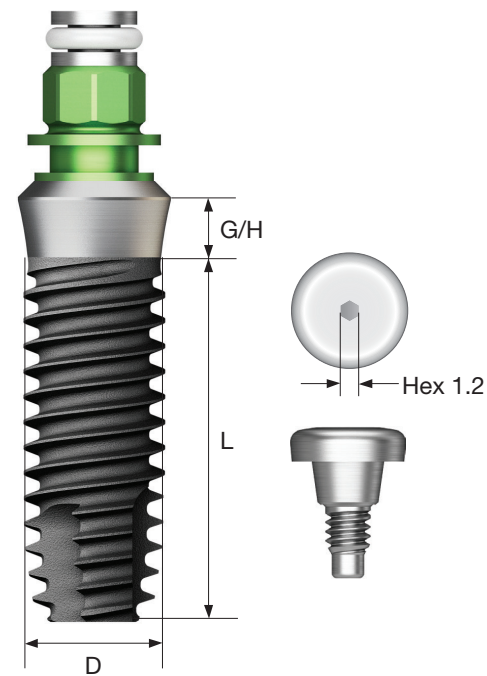

D Ø5.0
P Ø6.0
W

G/H \ L	6	7	8.5	10	11.5	13
	Short					
1.8	SS2W5006S18	SS2W5007S18	SS2W5008S18	SS2W5010S18	SS2W5011S18	SS2W5013S18
2.8	SS2W5006S28	SS2W5007S28	SS2W5008S28	SS2W5010S28	SS2W5011S28	SS2W5013S28

D Ø4.0
P Ø4.8
R

G/H \ L	7	8.5	10	11.5	13
1.8	SS2R4007S18	SS2R4008S18	SS2R4010S18	SS2R4011S18	SS2R4013S18
2.8		SS2R4008S28	SS2R4010S28	SS2R4011S28	SS2R4013S28

D Ø4.5
P Ø4.8
R

G/H \ L	7	8.5	10	11.5	13
1.8	SS2R4507S18	SS2R4508S18	SS2R4510S18	SS2R4511S18	SS2R4513S18
2.8		SS2R4508S28	SS2R4510S28	SS2R4511S28	SS2R4513S28

D Ø4.5
P Ø6.0
W

G/H \ L	7	8.5	10	11.5	13
1.8	SS2W4507S18	SS2W4508S18	SS2W4510S18	SS2W4511S18	SS2W4513S18
2.8	SS2W4507S28	SS2W4508S28	SS2W4510S28	SS2W4511S28	SS2W4513S28

Nominal diameter may differ from the actual diameter of the product
Note: Short implant should be used after a sufficient healing period. It is used by splinting with other implants for prosthesis

SSIII SA Implant

- Non-submerged type implant with an internal octa 8° tapered connection based on the one-stage surgery
- Optimal thread design for realization of optimal SA surface
- Tapered body design for excellent initial stability
- Superior self-threading effect with corkscrew thread
- Ensuring excellent primary stability needed for immediate loading even in soft bone

Ultra-wide


- Useful for placement in a fresh extraction socket in the posterior region, immediate placement case or for replacing a failed implant
 - Optimized apex design for excellent initial stability in a fresh extraction socket or in 3mm from the bottom
 - Recommended placement torque: $\leq 40\text{Ncm}$
- ※ Implants with a diameter of 4.5mm or greater are recommended for the posterior region with a single case

NoMount implant order code

: implant product code (ex : SS3R4011S18)

Pre-Mounted implant (implant + simple mount + cover screw) order code

: A + implant product code (ex : ASS3R4011S18)

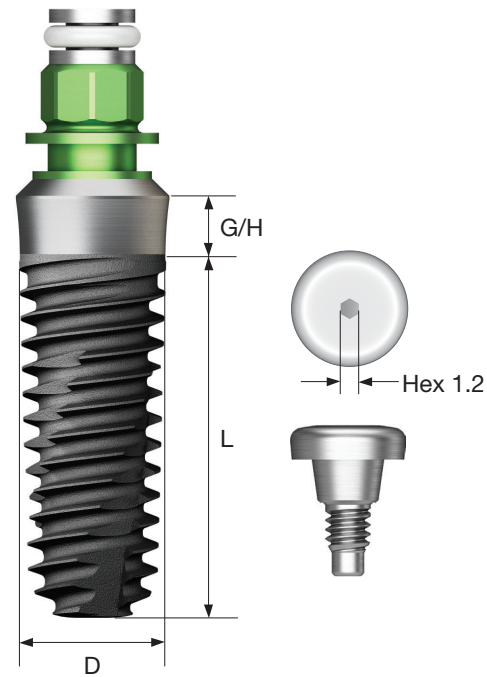

D Ø4.5
P Ø4.8
R

NEW 2025

G/H \ L	7	7	8.5	10	11.5	13
0.8	Extra short SS3R4506S08					
1.8		SS3R4507S18	SS3R4508S18	SS3R4510S18	SS3R4511S18	SS3R4513S18
2.8			SS3R4508S28	SS3R4510S28	SS3R4511S28	SS3R4513S28

D Ø4.5
P Ø6.0
W

G/H \ L	7	7	8.5	10	11.5	13
1.8	Extra short SS3W4506S18	SS3W4507S18	SS3W4508S18	SS3W4510S18	SS3W4511S18	SS3W4513S18
2.8		SS3W4507S28	SS3W4508S28	SS3W4510S28	SS3W4511S28	SS3W4513S28

D Ø3.5
P Ø4.8
R

G/H \ L	8.5	10	11.5	13
1.8	SS3R3508S18	SS3R3510S18	SS3R3511S18	SS3R3513S18
2.8	SS3R3508S28	SS3R3510S28	SS3R3511S28	SS3R3513S28

D Ø4.0
P Ø4.8
R

NEW 2025

G/H \ L	7	7	8.5	10	11.5	13
0.8	Extra short SS3R4006S08					
1.8	SS3R4006S18	SS3R4007S18	SS3R4008S18	SS3R4010S18	SS3R4011S18	SS3R4013S18
2.8			SS3R4008S28	SS3R4010S28	SS3R4011S28	SS3R4013S28

D Ø5.0
P Ø6.0
W

NEW 2025

G/H \ L	6	6	6	7	8.5	10	11.5	13
0.8	Extra short SS3W5004S08	Extra short SS3W5005S08	Short					
1.8		SS3W5005S18	SS3W5006S18	SS3W5007S18	SS3W5008S18	SS3W5010S18	SS3W5011S18	SS3W5013S18
2.8			SS3W5006S28	SS3W5007S28	SS3W5008S28	SS3W5010S28	SS3W5011S28	SS3W5013S28







Nominal diameter may differ from the actual diameter of the product

Note: Short implant should be used after a sufficient healing period. It is used by splinting with other implants for prosthesis







SSIII SA Implant

Ultra-wide

D \varnothing 6.0
P \varnothing 6.0
W

G/H \ L	6	7	8.5	10	11.5	13
	 Short					
1.8	SS3W6006S18	SS3W6007S18	SS3W6008S18	SS3W6010S18	SS3W6011S18	SS3W6013S18
2.8	SS3W6006S28	SS3W6007S28	SS3W6008S28	SS3W6010S28	SS3W6011S28	SS3W6013S28

D \varnothing 7.0
P \varnothing 6.0
W

G/H \ L	6	7	8.5	10	11.5	13
	 Short					
1.8	SS3W7006S18	SS3W7007S18	SS3W7008S18	SS3W7010S18	SS3W7011S18	SS3W7013S18
2.8	SS3W7006S28	SS3W7007S28	SS3W7008S28	SS3W7010S28	SS3W7011S28	SS3W7013S28

OSSTEM[®]
IMPLANT

Nominal diameter may differ from the actual diameter of the product

Note Short implant should be used after a sufficient healing period. It is used by splinting with other implants for prosthesis

Mount & Screw

Simple Mount NEW 2025

- Selected according to the implant platform
- Hand tightened with a 1.2 hex driver
- ※ Disposable, Do not reuse
- P = Platform

R Regular

W Wide

P

R



SSHRG

W



SSHWB

Closing Screw NEW 2025

- For lack of soft tissue in the suture
- Hand tightened with a 1.2 hex driver
- P = Platform

R Regular

W Wide

P

R



SSCS480N

W



SSCS600N

Cover Screw


- Selected according to the implant platform
- Hand tightened with a 1.2 hex driver
- P = Platform

R Regular

W Wide

P

R

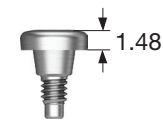

SSCS480

W




SSCS600

Healing Abutment

- Selected according to the implant platform
- Hand tightened with a 1.2 hex driver
- P = Platform

- R** Regular
- W** Wide

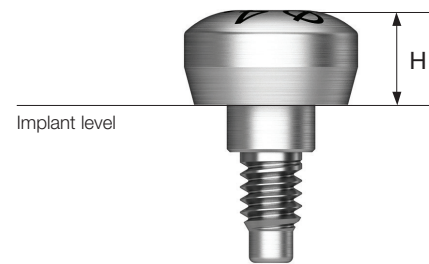

OSSTEM[®]
IMPLANT

098

P \ H	2.0	3.0	4.0	5.0
Ø4.8	SSH482	SSH483	SSH484	SSH485

P \ H	2.0	3.0	4.0	5.0
Ø6.0	-	SSH603	SSH604	SSH605


099

SS SYSTEM

SS SYSTEM

Solid / Excellent Solid


Abutment Level Impression

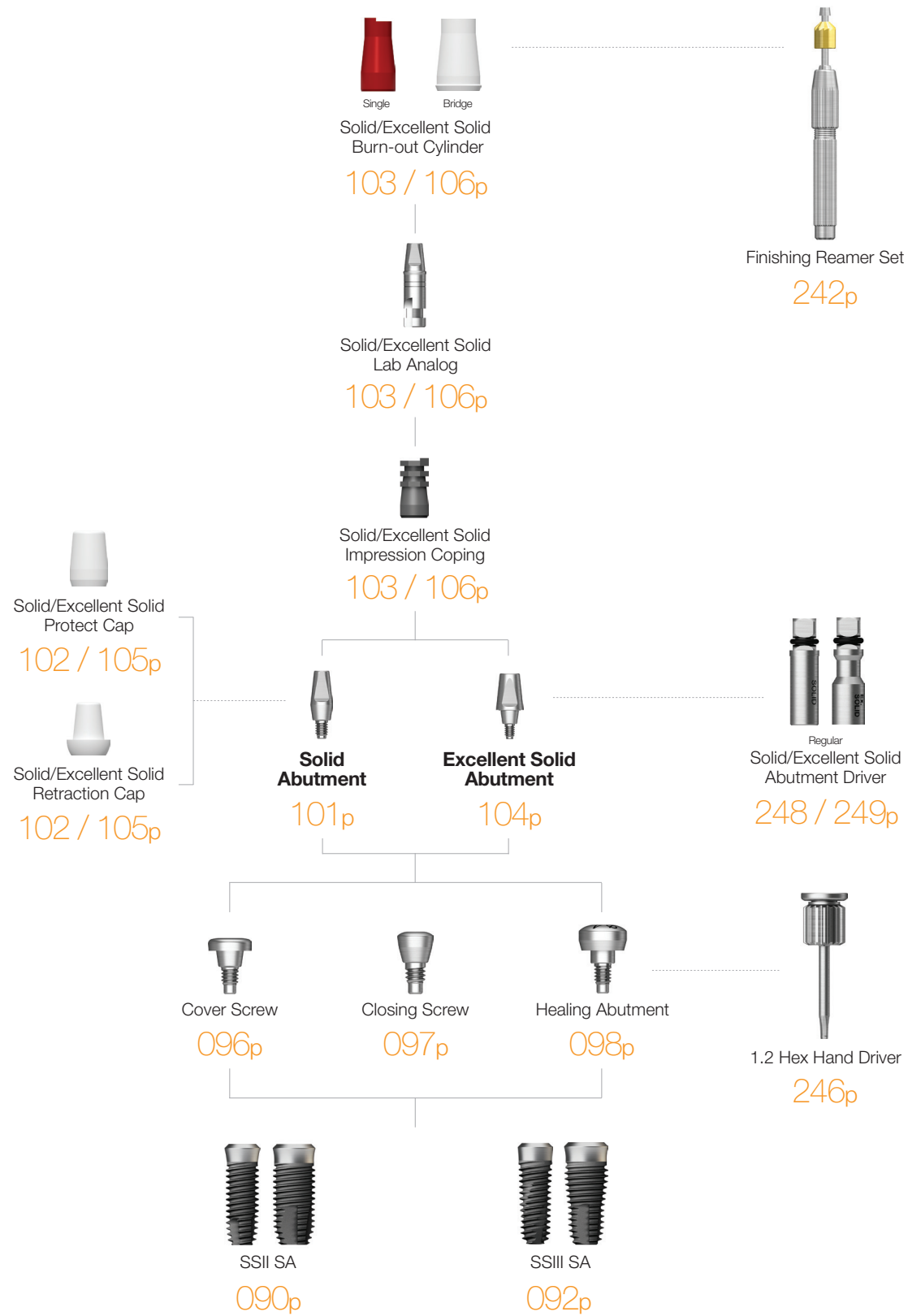

Solid Abutment

- Abutment for producing cement-retained prosthesis
- Abutment level impression
- Ø4.8 : Tightened with a Solid Abutment driver (code : SDSL/SDSS)
- Ø6.0 : Tightened with a 1.2 hex driver or Solid Abutment driver (code : SD60S)
- Recommended tightening torque : 30Ncm
- Packing unit : abutment + protect cap


Abutment + protect cap order code
: product code + **P** (ex : SSS485**P**)

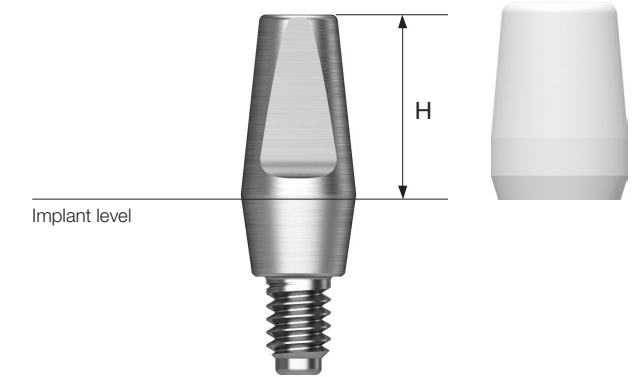

P Ø4.8



P Ø6.0

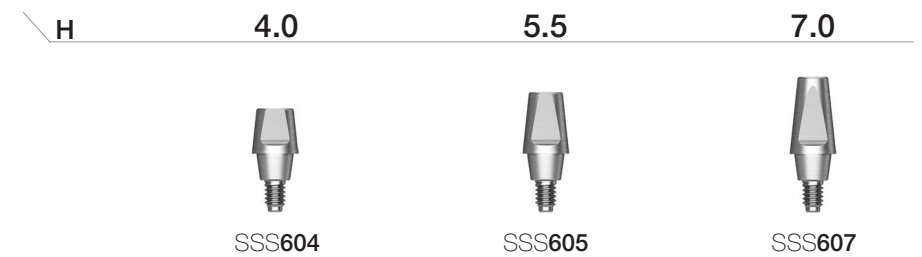







Solid Abutment Components

Solid Protect Cap

- Solid Abutment protection with reduced patient discomfort
- Used as a temporary crown base

R Regular

W Wide







P \ H	4.0	5.5	7.0
Ø 4.8	 SSC484	 SSC485	 SSC487
Ø 6.0	 SSC604	 SSC605	 SSC607

Solid Impression Coping

- Components for Solid Abutment impression
- Enabling production of elaborate prosthesis using lab analog
- Used by selecting the color matching the abutment height

R Regular

W Wide







P \ H	4.0	5.5	7.0
Ø 4.8	 SSIC484	 SSIC485	 SSIC487
Ø 6.0	 SSIC604	 SSIC605	 SSIC607

Solid Retraction Cap

- Ensuring clear margin by pushing the gingiva around the margin in the direct impression of Solid Abutment
- Used as a temporary crown base

R Regular

W Wide







P \ H	4.0	5.5	7.0
Ø 4.8	 SSSRC484	 SSSRC485	 SSSRC487
Ø 6.0	 SSSRC604	 SSSRC605	 SSSRC607

Solid Lab Analog

- Components for Solid Abutment reproduction on a model after impression taking
- Used by assembling to the solid impression coping in the same color

R Regular

W Wide

P \ H	4.0	5.5	7.0
Ø 4.8	 SSSA484	 SSSA485	 SSSA487
Ø 6.0	 SSSA604	 SSSA605	 SSSA607

Solid Burn-out Cylinder

- Components replacing the resin cap prior to fabrication of wax-up using Solid Abutment
- Enabling the production of elaborate prosthesis with uniform interior
- Used after removing the tightening connection of lower margin after casting

R Regular


W Wide

P \ H	Single	Bridge
Ø 4.8	 SSSP480S	 SSSP480B
Ø 6.0	 SSSP600S	 SSSP600B

Excellent Solid Abutment

- Abutment for producing cement-retained prosthesis
- Larger in volume compared to Solid Abutment, suitable for molars or when removal is required
- Abutment level impression
- Ø 4.8 : Tightened with a 1.2 hex driver or Excellent Solid Abutment driver (code : ESDSS/ESDSL)
- Ø 6.0 : Tightened with a 1.2 hex driver or Excellent Solid Abutment driver (code : ESD60S)
- Recommended tightening torque: 30Ncm
- Packing unit : abutment + protect cap

Abutment + protect cap order code
: product code + **P** (ex : SSE485**P**)

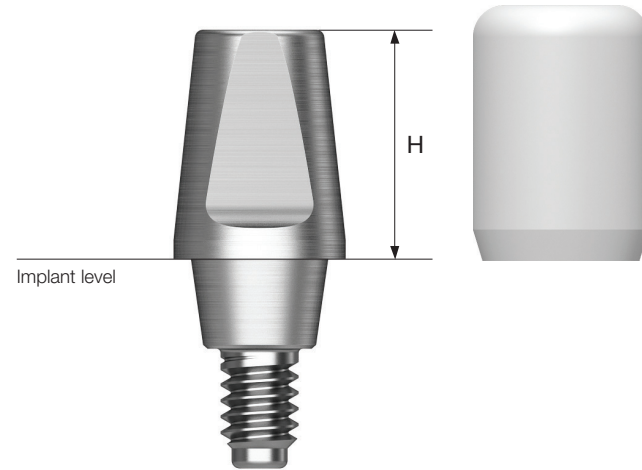

Excellent Solid Abutment Components

Excellent Solid Protect Cap

- Used for Excellent Solid Abutment protection and reducing patient discomfort
- Used as a temporary crown base

R Regular
W Wide

P \ H	4.0	5.5	7.0
Ø 4.8	SSEC484	SSEC485	SSEC487
Ø 6.0	SSEC604	SSEC605	SSEC607

Excellent Solid Retraction Cap

- Used for accurate margin reproduction by pushing away the surrounding gingiva when taking a direct impression of Excellent Solid Abutment
- Used as a temporary crown base

R Regular
W Wide

P \ H	4.0	5.5	7.0
Ø 4.8	SSERC484	SSERC485	SSERC487
Ø 6.0	SSERC604	SSERC605	SSERC607

P Ø 4.8

R

H	4.0	5.5	7.0
	SSE484	SSE485	SSE487

P Ø 6.0

W




H	4.0	5.5	7.0
	SSE604	SSE605	SSE607

Excellent Solid Abutment Components

Excellent Solid Impression Coping

- Components for Excellent Solid Abutment impression
- Enabling production of elaborate prosthesis using lab analog
- Used by selecting the color matching the abutment height

- R** Regular
- W** Wide

P \ H	4.0	5.5	7.0
			
	SSEIC484	SSEIC485	SSEIC487
	SSEIC604	SSEIC605	SSEIC607

Ø 4.8
Ø 6.0

Excellent Solid Lab Analog

- Components for Excellent Solid Abutment reproduction on model after impression taking
- Used by connecting to the appropriate color coded Excellent Solid Impression Cap

- R** Regular
- W** Wide

P \ H	4.0	5.5	7.0
			
	SSEA484	SSEA485	SSEA487
	SSEA604	SSEA605	SSEA607

Ø 4.8
Ø 6.0

Excellent Solid Burn-out Cylinder

- Components replacing the resin cap prior to fabrication of wax-up using Excellent Solid Abutment
- Enabling the production of elaborate prosthesis with uniform interior
- Used after removing the tightening connection of lower margin after casting

- R** Regular
- W** Wide


P \ H	Single	Bridge
		
	SSEP480S	SSEP480B
	SSEP600S	SSEP600B

Ø 4.8
Ø 6.0

OSSTEM[®]
IMPLANT


ComOcta

Implant Level Impression

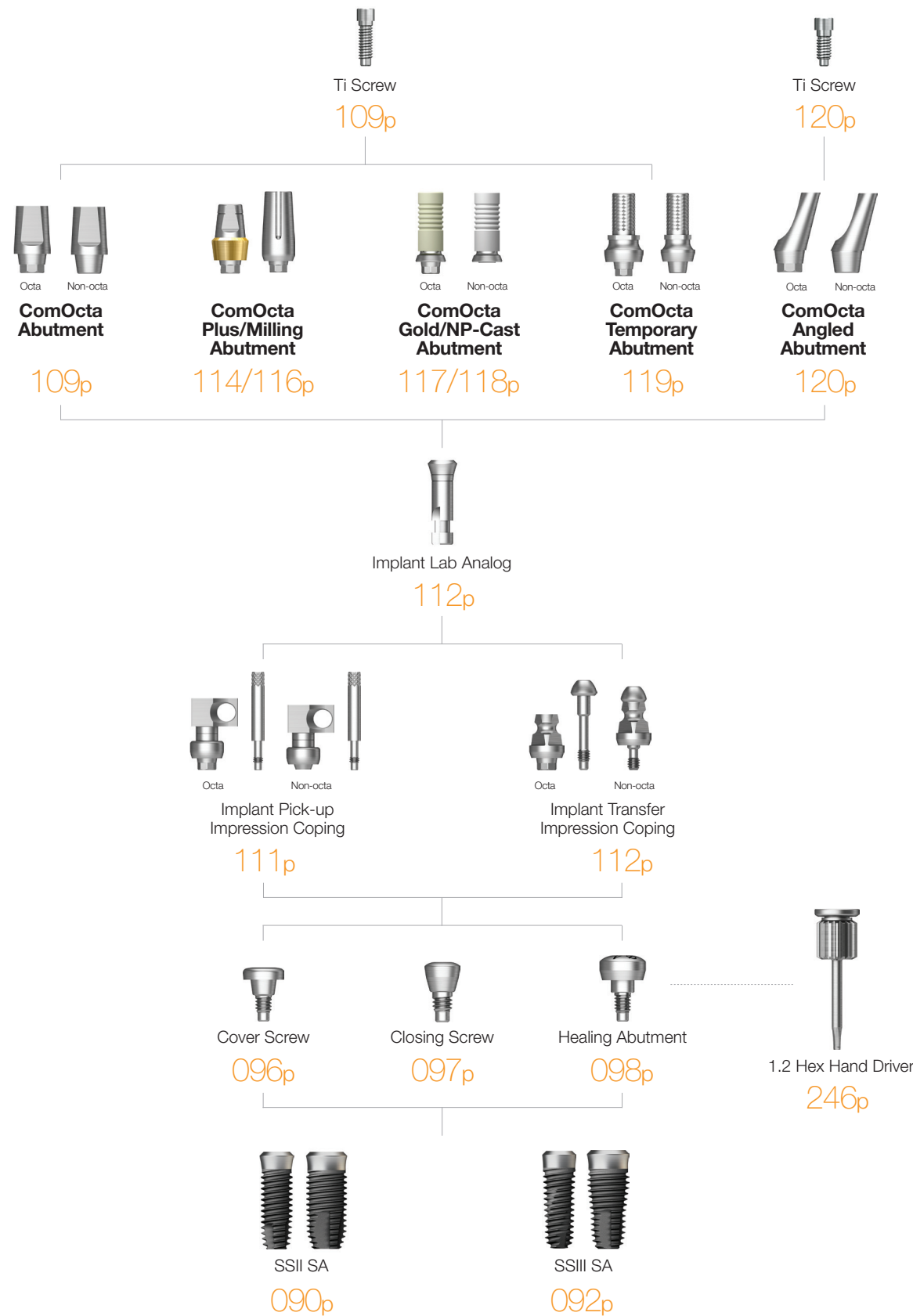

ComOcta Abutment

- Abutment for producing cement-retained / combination prosthesis
- Implant level impression
- Enabling abutment level impression using the retraction cap
- Tightened with a 1.2 hex driver
- Recommended tightening torque: 30Ncm
- Packing unit : abutment + Ti screw

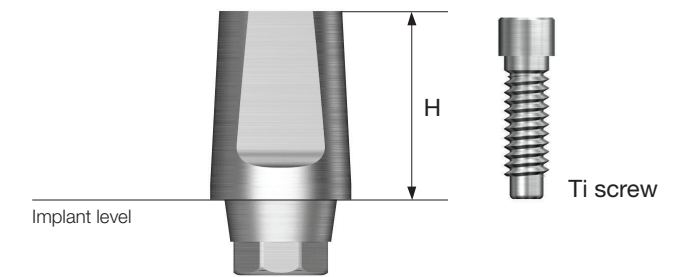

Abutment + Ti screw order code
: product code + TH (ex : SSCA485TH)

P Ø4.8



Ti screw : ASR200

H	4.0	5.5	7.0	4.0	5.5	7.0
Type	Octa			Non-Octa		
	SSCA484	SSCA485	SSCA487	SSCA484N	SSCA485N	SSCA487N

P Ø6.0



Ti screw : ASR200







H	4.0	5.5	7.0	4.0	5.5	7.0
Type	Octa			Non-Octa		
	SSCA604	SSCA605	SSCA607	SSCA604N	SSCA605N	SSCA607N

ComOcta Abutment Components

ComOcta Protect Cap NEW 2025

- ComOcta Abutment protection with reduced patient discomfort
- Used as a temporary crown base
- Excellent Solid Protect Cap used in common for wide type







R Regular
W Wide

P \ H	4.0	5.5	7.0
Ø 4.8	 SSCC484	 SSCC485	 SSCC487
Ø 6.0	 SSEC604	 SSEC605	 SSEC607

ComOcta Retraction Cap

- Used for accurate margin reproduction by pushing away the surrounding gingiva when taking a direct impression of ComOcta Abutment
- Used as a temporary crown base







R Regular
W Wide

P \ H	4.0	5.5	7.0
Ø 4.8	 SSCRC484	 SSCRC485	 SSCRC487
Ø 6.0	 SSCRC604	 SSCRC605	 SSCRC607

ComOcta Impression Coping

- Components for ComOcta Abutment impression
- Enabling production of elaborate prosthesis using lab analog
- Used by selecting the color matching the abutment height
- Excellent Solid Impression Coping used in common for wide type

R Regular
W Wide

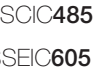P \ H	4.0	5.5	7.0
Ø 4.8	 SSCIC484	 SSCIC485	 SSCIC487
Ø 6.0	 SSEIC604	 SSEIC605	 SSEIC607

ComOcta Lab Analog

- Components for Excellent Solid Abutment reproduction on model after impression taking
- Used by connecting to the same color as the ComOcta Impression Cap


R Regular
W Wide










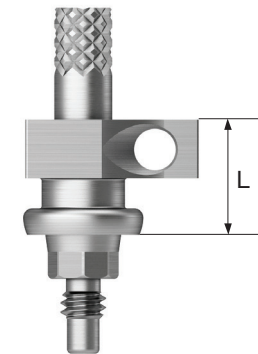

P \ H	4.0	5.5	7.0
Ø 4.8	 SSCLA484	 SSCLA485	 SSCLA487
Ø 6.0	 SSCLA604	 SSCLA605	 SSCLA607

Implant Pick-up Impression Coping

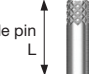- Components for implant level impression taking
- Using open tray
- Unique design stably fixed within the impression body
- Hand tightened with a 1.2 hex driver
- Packing unit : impression coping body + guide pin(*)

R Regular (Green)
W Wide (Blue)




P \ L	Type	5		10		Guide Pin		
		Octa	Non-Octa	Octa	Non-Octa	10	15	17
Ø 4.8								
Ø 6.0								
		SSICAS480	SSICAS480N	SSICA480	SSICA480N	CSR100*(L5) CSR150*(L10) CSR170		
		SSICAS600	SSICAS600N	SSICA600	SSICA600N			

ComOcta Abutment Components

Implant Transfer Impression Coping

- Components for implant level impression taking
- Using closed tray
- Triangular arc structure for stable fastening and accurate repositioning
- Hand tightened with a 1.2 hex driver
- Packing unit
 - Octa : impression coping body + guide pin
 - Non-octa : impression coping

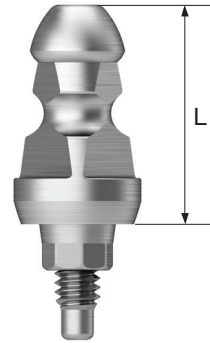

- R** Regular (Green)
- W** Wide (Blue)

P \ L	9.5		12.5		9.5		12.5	
	Type				Non-Octa			
	Octa				Non-Octa			
	Ø 4.8	Ø 6.0	SSCTIS480	SSCTIS600	SSCTIL480	SSCTIL600	SSCTIS480N	SSCTIS600N
			SSCTIL480	SSCTIL600	SSCTIL480N	SSCTIL600N		



Implant Lab Analog

- Lab analog for implant level impression
- Used by selecting according to the implant platform Ø4.8/6.0

- R** Regular (Green)
- W** Wide (Blue)

P \ L	Type
	SSFA480
	SSFA600

Ø 4.8
Ø 6.0




SSFA480
SSFA600

ComOcta Plus Abutment

- Abutment for producing cement-retained/combination prosthesis
- Used for thick gingiva or deeply placed implant
- 45° platform contact for abutment-implant connection
- Implant level impression
- Tightened with a 1.2 hex driver
- Recommended tightening torque: 30Ncm
- Packing unit : abutment + Ti screw

Abutment + Ti screw order code
: product code + **TH** (ex : SSCAP4826**CTH**)

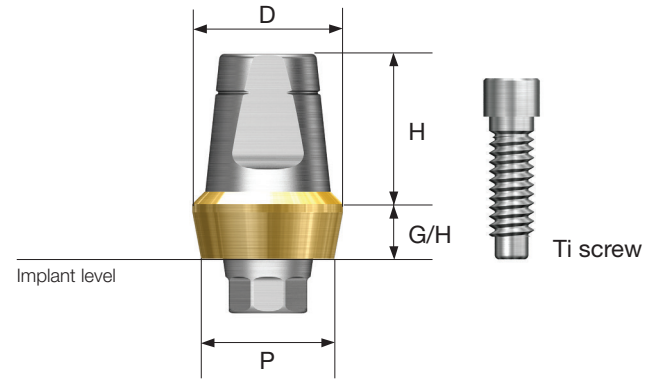

P Ø6.0



Ti screw
: ASR200

		H \ G/H D	1.0 Ø6.5	2.0 Ø6.5	3.0 Ø6.5	4.0 Ø6.5
Octa	4.0		SSCAP6014C	SSCAP6024C	SSCAP6034C	SSCAP6044C
	5.5		SSCAP6016C	SSCAP6026C	SSCAP6036C	SSCAP6046C
Non-Octa	5.5		SSCAP6016CN	SSCAP6026CN	SSCAP6036CN	SSCAP6046CN

		H \ G/H D	2.0 Ø6.8	3.0 Ø7.2	4.0 Ø7.6
--	--	--------------	-------------	-------------	-------------

Octa	5.5		SSCAP6026E	SSCAP6036E	SSCAP6046E
Non-Octa	5.5		SSCAP6026EN	SSCAP6036EN	SSCAP6046EN

114

P Ø4.8



Ti screw
: ASR200

		H \ G/H D	1.0 Ø5.5	2.0 Ø5.5	3.0 Ø5.5	4.0 Ø5.5
Octa	4.0		SSCAP4814C	SSCAP4824C	SSCAP4834C	SSCAP4844C
	5.5		SSCAP4816C	SSCAP4826C	SSCAP4836C	SSCAP4846C
Non-Octa	5.5		SSCAP4816CN	SSCAP4826CN	SSCAP4836CN	SSCAP4846CN

		H \ G/H D	2.0 Ø6.0	3.0 Ø6.5	4.0 Ø7.0
--	--	--------------	-------------	-------------	-------------

Octa	5.5		SSCAP4826E	SSCAP4836E	SSCAP4846E
Non-Octa	5.5		SSCAP4826EN	SSCAP4836EN	SSCAP4846EN


115

SS SYSTEM

ComOcta Milling Abutment

- Abutment for producing cement-retained / combination prosthesis
- Used for customizing the shape of the abutment margin
- 45° platform contact for abutment-implant connection
- Implant level impression
- Tightened with a 1.2 hex driver
- Recommended tightening torque: 30Ncm
- Packing unit : abutment + Ti screw


Abutment + Ti screw order code
: product code + **TH** (ex : SSCMA4830**TH**)

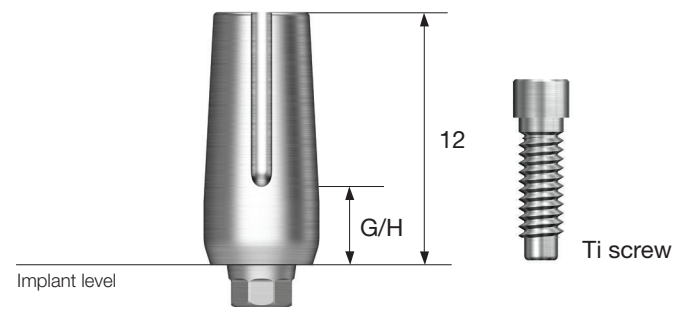







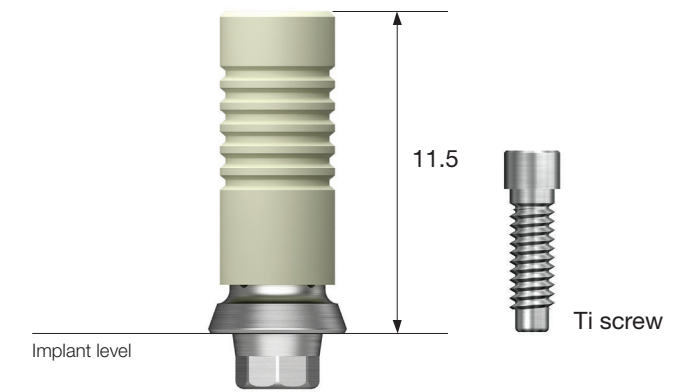







ComOcta Gold Abutment

- Abutment for producing cement-retained / combination / screw-retained prosthesis
- Used for fabrication of customized prosthesis by casting with gold alloys
- 45° platform contact for abutment-implant connection
- Abutment melting temperature: 1,400~1,450°C
- Implant level impression
- Tightened with a 1.2 hex driver
- Recommended tightening torque: 30Ncm
- Packing unit : abutment + Ti screw

Abutment + Ti screw order code
: product code + **TH** (ex : COG480**STH**)




<p>P Ø4.8 / G/H</p> <p>R</p> <p>Ti screw : ASR200</p>	<p>3.0</p>  <p>SSCMA4830</p>	<p>P Ø6.0 / G/H</p> <p>W</p> <p>Ti screw : ASR200</p>	<p>3.0</p>  <p>SSCMA6030</p>
---	---	---	---

<p>P Ø4.8 / Type</p> <p>R</p> <p>Ti screw : ASR200</p>	<table border="0"> <tr> <td>Octa</td> <td>Non-Octa</td> </tr> <tr> <td></td> <td></td> </tr> <tr> <td>COG480S</td> <td>COG480B</td> </tr> </table>	Octa	Non-Octa			COG480S	COG480B	<p>P Ø6.0 / Type</p> <p>W</p> <p>Ti screw : ASR200</p>	<table border="0"> <tr> <td>Octa</td> <td>Non-Octa</td> </tr> <tr> <td></td> <td></td> </tr> <tr> <td>COG600S</td> <td>COG600B</td> </tr> </table>	Octa	Non-Octa			COG600S	COG600B
Octa	Non-Octa														
															
COG480S	COG480B														
Octa	Non-Octa														
															
COG600S	COG600B														

ComOcta NP-Cast Abutment

- Abutment for producing cement-retained/combination/screw-retained prosthesis
- Used for fabrication of customized prosthesis by casting with nonprecious metal alloys
- 45° platform contact for abutment-implant connection
- Abutment melting temperature: 1,400~1,450°C
- Implant level impression
- Tightened with a 1.2 hex driver
- Recommended tightening torque: 30Ncm
- Packing unit : abutment + Ti screw


Abutment + Ti screw order code
: product code + **TH** (ex : CON480**STH**)

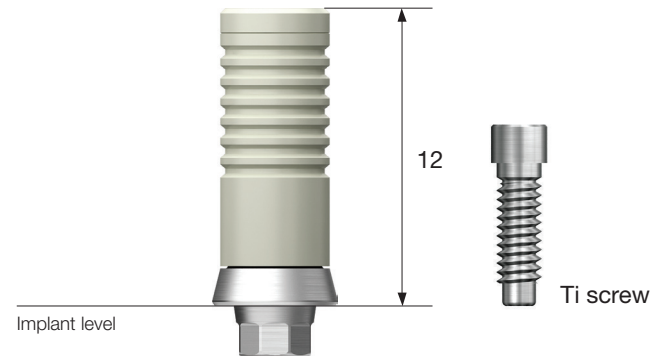

ComOcta Temporary Abutment

- Abutment for producing cement-retained/screw-retained temporary prosthesis
- Used by removing for producing temporary prosthesis (Ti Gr-3)
- Implant level impression
- Tightened with a 1.2 hex driver
- Recommended tightening torque: 20Ncm
- Packing unit : abutment + Ti screw

Abutment + Ti screw order code
: product code + **TH** (ex : SSTAO480**TH**)

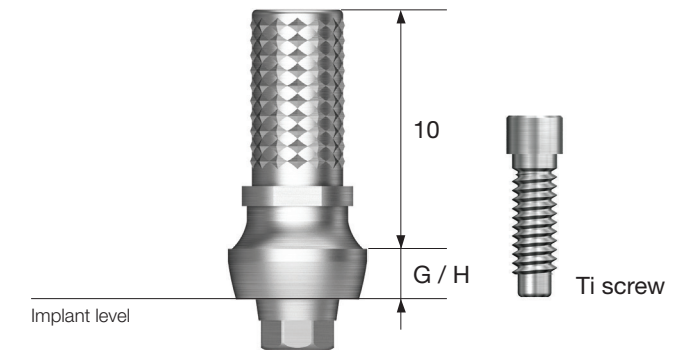

P Ø4.8 Type **Octa** **Non-Octa**

R
Ti screw : ASR200

CON480S CON480B

P Ø6.0 Type **Octa** **Non-Octa**

W
Ti screw : ASR200

CON600S CON600B

P Ø4.8 G/H Type **0** **2.0** **0** **2.0**

R **Octa** **Non-Octa**
Ti screw : ASR200

SSTAO480 SSTAO482 SSTAN480 SSTAN482


P Ø6.0 G/H Type **0** **2.0** **0** **2.0**

W **Octa** **Non-Octa**
Ti screw : ASR200

SSTAO600 SSTAO602 SSTAN600 SSTAN602

ComOcta Angled Abutment

- Abutment for producing cement-retained/combination prosthesis
- 15°/25° implant placement angle compensation
- Dedicated abutment screws are used
- Implant level impression
- Tightened with a 1.2 hex driver
- Recommended tightening torque: 30Ncm
- Packing unit : abutment + Ti screw(only angled)

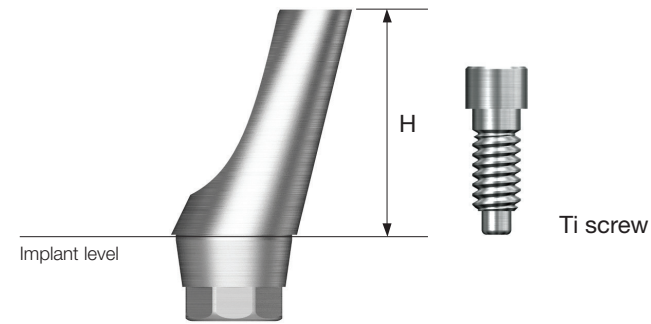

Abutment + Ti screw order code
: product code + **TH** (ex : SSA4815**TH**)

OSSTEM[®]
IMPLANT

P Ø4.8



Ti screw
: ASS200

H	Angle Type	15° Octa	20°	15° Non-Octa	20°
6.7					
		SSA4815	SSA4820	SSA4815N	SSA4820N

P Ø6.0




Ti screw
: ASS200

H	Angle Type	15° Octa	20°	15° Non-Octa	20°
5.5					
		SSA6015	SSA6020	SSA6015N	SSA6020N


Octa

Abutment Level Impression

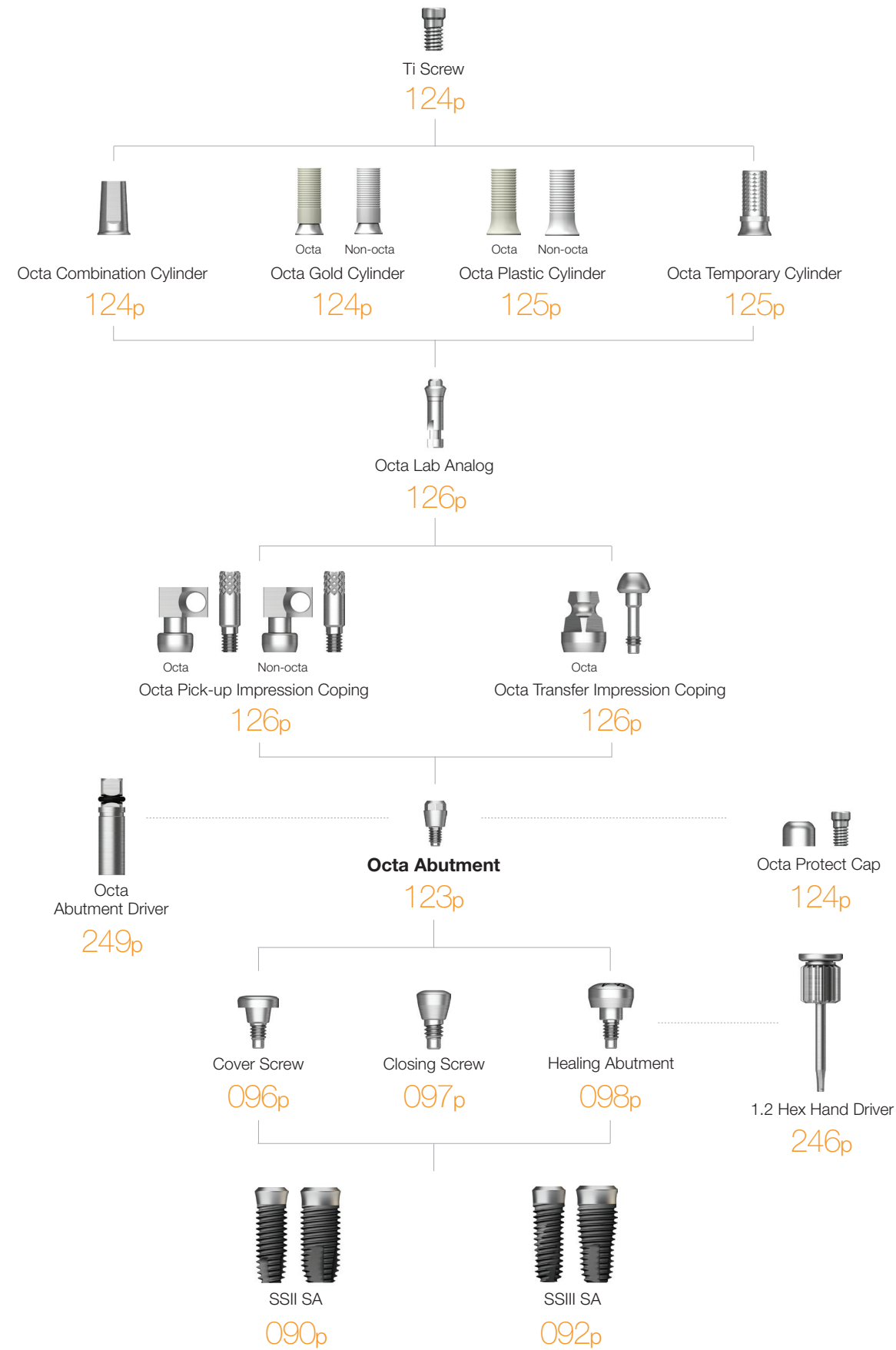

Octa Abutment

- Used for producing combination/screw-retained prosthesis in multiple case
- Implant placement angle compensated up to 60°
- Tightened with a dedicated outer driver (code : ODSL/ODSS)
- Recommended tightening torque: 30Ncm

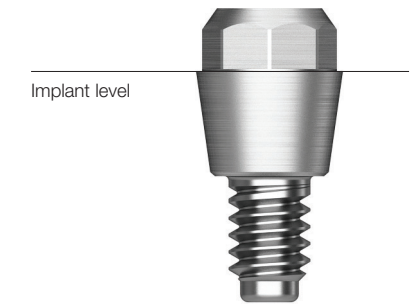

P Ø4.8



SSOA480

P Ø6.0



SSOA600

Octa Abutment Components

Octa Protect Cap

- Protect cap for Octa Abutment
- Hand tightened with a 1.2 hex driver
- Packing unit : protect cap + Ti screw

Protect cap + Ti screw order code
: product code + **TH** (ex : SSHC480**TH**)

R Regular

W Wide




Octa Temporary Cylinder

- Used for producing temporary prosthesis in Octa Abutment (Ti Gr-3)
- Tightened with a 1.2 hex driver
- Recommended tightening torque: 20Ncm
- Packing unit : cylinder + Ti cylinder screw

Cylinder + Ti screw order code
: product code + **TH** (ex : SSTCO480**TH**)

R Regular

W Wide

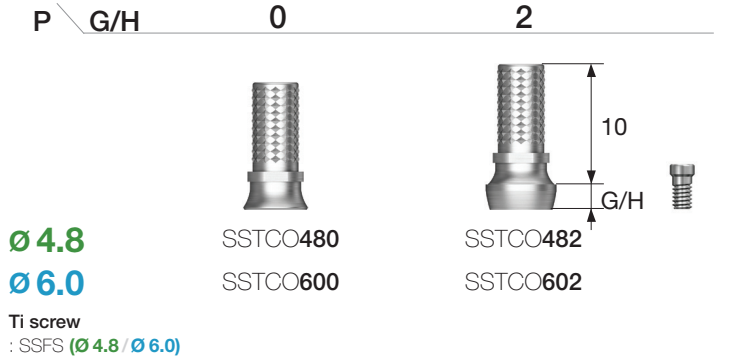


Octa Gold Cylinder

- Used for producing screw-retained prosthesis in Octa Abutment
- Used for fabrication of customized prosthesis by casting with gold alloys
- Cylinder melting temperature: 1,400~1,450°C
- Tightened with a 1.2 hex driver
- Recommended tightening torque: 20Ncm
- Packing unit : cylinder + Ti cylinder screw

Cylinder + Ti screw order code
: product code + **TH** (ex : SSGCO480**TH**)

R Regular

W Wide

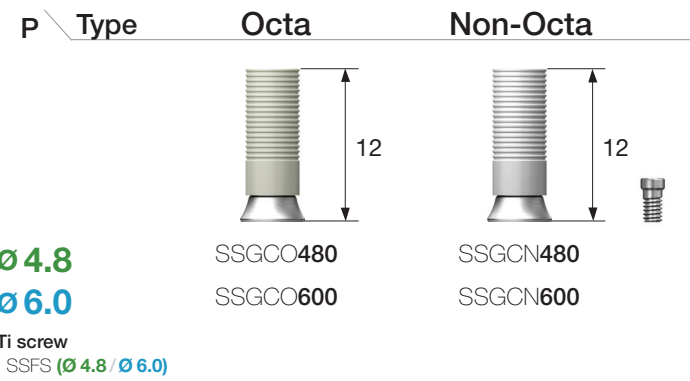


Octa Plastic Cylinder

- Used for producing screw-retained prosthesis in Octa Abutment
- Used for fabrication of customized prosthesis by casting with nonprecious metal alloys
- Tightened with a 1.2 hex driver
- Recommended tightening torque: 20Ncm
- Packing unit : cylinder + Ti cylinder screw

Cylinder + Ti screw order code
: product code + **TH** (ex : SSPSO480**TH**)

R Regular

W Wide

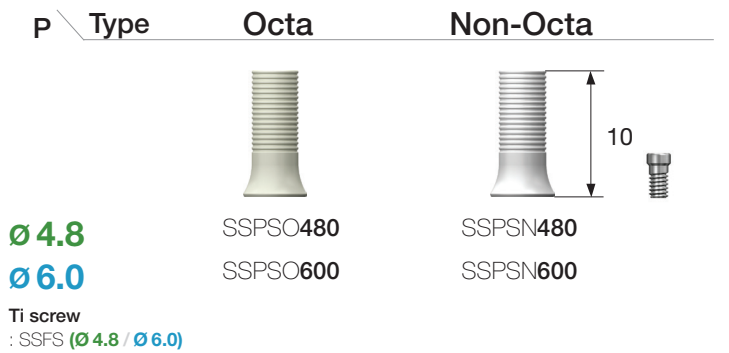


Octa Combination Cylinder

- Used for producing combination prosthesis in Octa Abutment
- Connection structure for both octa/non-octa
- Tightened with a 1.2 hex driver
- Recommended tightening torque: 20Ncm
- Packing unit : cylinder + Ti cylinder screw

Cylinder + Ti screw order code
: product code + **TH** (ex : SSOCC480**TH**)

R Regular

W Wide

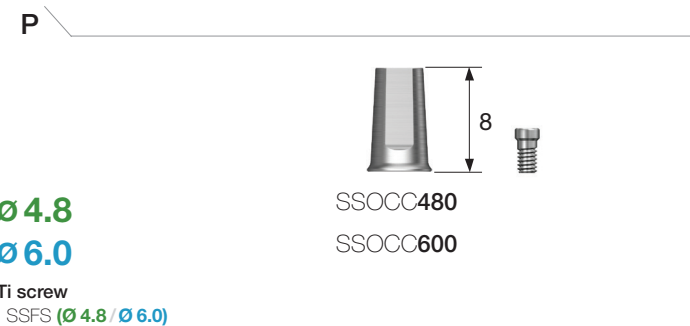


Octa Abutment Components

Octa Pick-up Impression Coping

- Pick up impression coping for Octa Abutment
- Hand tightened with a 1.2 hex driver
- Packing unit : impression coping body + guide pin(*)

R Regular (Green)

W Wide (Blue)

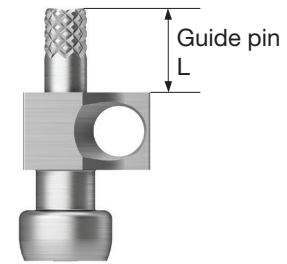

P \ L	Octa		Non-Octa		Guide Pin	
	Ø 4.8	Ø 6.0	Ø 4.8	Ø 6.0	0	5.0
	SSICO480	SSICO600	SSICN480	SSICN600	SSGS100	SSGS150*

OSSTEM[®]
IMPLANT

Octa Transfer Impression Coping

- Transfer impression coping for Octa Abutment
- Hand tightened with a 1.2 hex driver
- Packing unit : Impression coping body + guide pin

R Regular (Green)

W Wide (Blue)

P	Ø 4.8	Ø 6.0
	SSOTI480	SSOTI600

Octa Lab Analog

- Lab analog for Octa Abutment
- Hand tightened with a 1.2 hex driver


R Regular (Green)

W Wide (Blue)

P	Ø 4.8	Ø 6.0
	SSLA480	SSLA600


Port / Locator® / O-ring

Overdenture

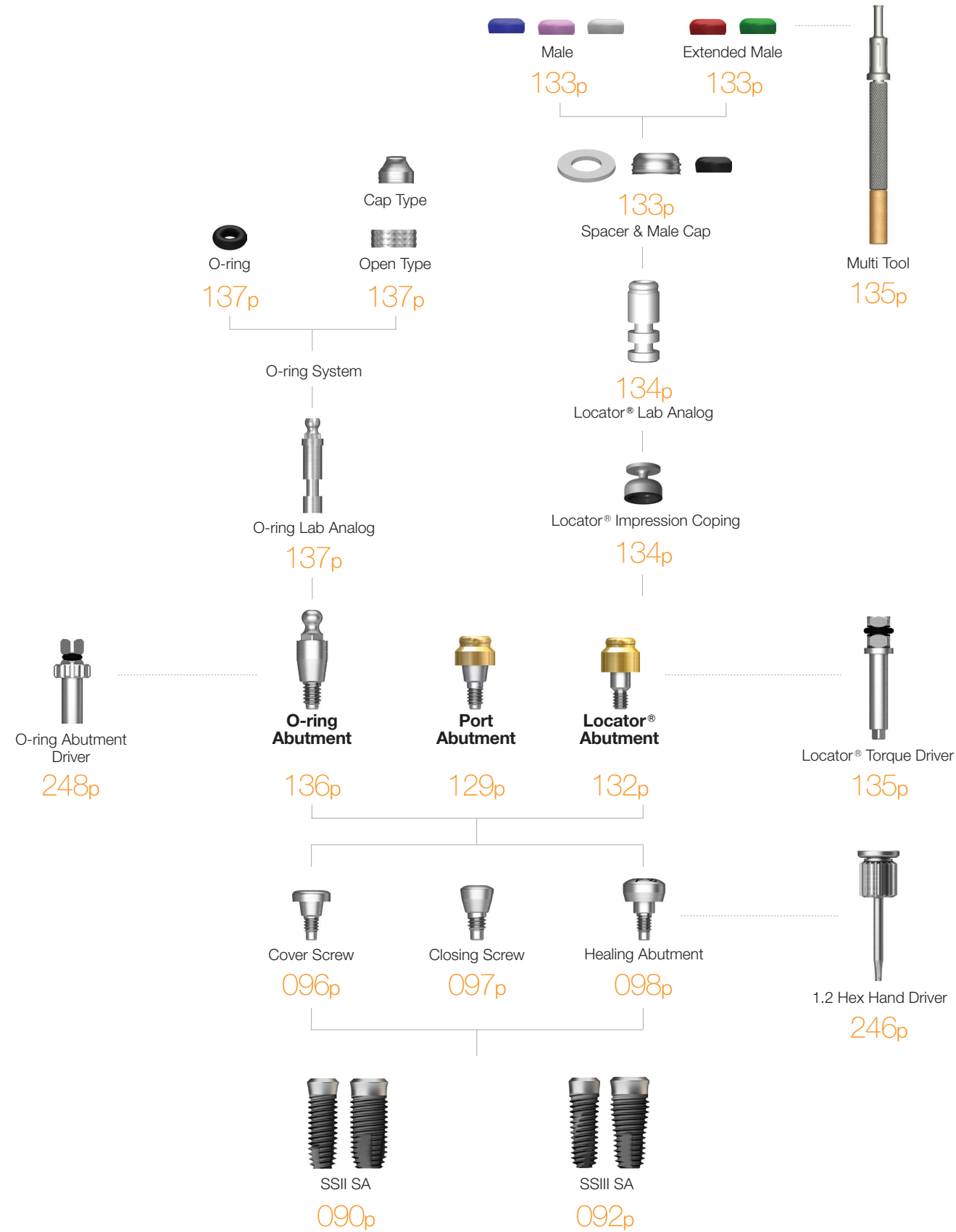


Port Abutment

- Placement angle compensated up to 40°
- Vertical dimension lower by 1.5mm, construction of various attachments with stable fixing
- Tightened with a dedicated outer driver (code : TWLDLK/TWLDLSK)
- Recommended tightening torque: 30Ncm

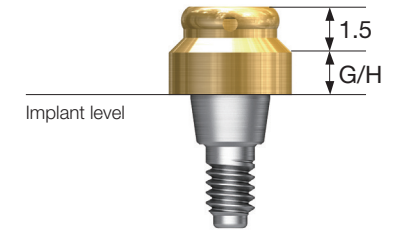

P Ø4.8

R

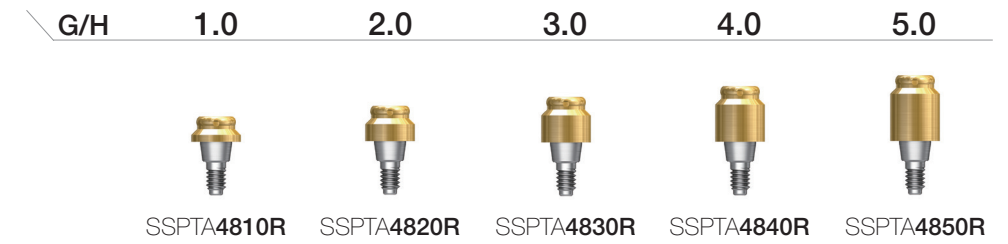

Port Abutment Components

Port Male KIT

- Components
 - Block out spacer / denture cap connected black processing male
 - Replacement male blue/pink/clear
- Used by selecting the male with adequate retention force for each case
- Using a locator core tool for replacing the male
- Packing unit : 1set



Port Replacement Male

- Retention force : Low retention (Approx. 6N)
- Placement angle compensated up to 20°
- Packing unit : 4ea



- Retention force : Medium retention (Approx. 12N)
- Placement angle compensated up to 20°
- Packing unit : 4ea



- Retention force : High retention (Approx. 22N)
- Placement angle compensated up to 20°
- Packing unit : 4ea



Port Extended Replacement Male

- Retention force : Low retention (Approx. 6N)
- Placement angle compensated up to 20~40°
- Packing unit : 4ea



- Retention force : Medium retention (Approx. 12N)
- Placement angle compensated up to 20~40°
- Packing unit : 4ea



Port Black Processing Male

- Male used only in prosthesis fabrication process
- Packing unit : 4ea



Port Male Cap

- Fixed to the denture by connecting with the male
- Packing unit : 1ea



Port Block Out Spacers

- Used for sealing of the space between the abutment and the denture cap when attaching the overdenture and denture cap in the oral cavity
- Packing unit : 20ea




Port Impression Coping

- Pick up impression coping for Locator Abutment
- Using closed tray
- Packing unit : Impression coping + Provisional male 1set



Port Lab Analog

- Lab analog for Locator Abutment
- Packing unit : 1ea

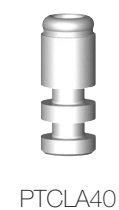


Port Multi Tool

- Used for placing and removing the replacement male in the denture cap
- Separated into three pieces and used as a hand driver for port abutment



Locator® Abutment

- Locator Abutment for SS abutment of Zest Dental Solutions
- Placement angle compensated up to 40°
- Vertical dimension lower by 1.5mm, construction of various attachments with stable fixing
- Tightened with a dedicated outer driver (code : TWLDLK/TWLDLSK)
- Recommended tightening torque: 30Ncm

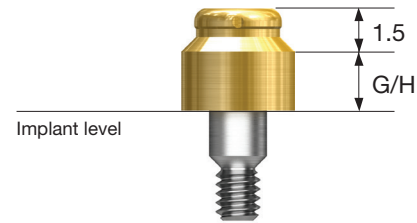

P Ø4.8



Locator® Abutment Components

Locator® Male Processing Kit

- Components
 - Block out spacer / denture cap connected black processing male
 - Replacement male blue/pink/clear
- Used by selecting the male with adequate retention force for each case
- Using a locator core tool for replacing the male
- Packing unit : 2set



Locator® Replacement Male

- Retention force : Approx. 6N
- Placement angle compensated up to 20°
- Packing unit : 4ea

- Retention force : Approx. 12N
- Placement angle compensated up to 20°
- Packing unit : 4ea

- Retention force : Approx. 22N
- Placement angle compensated up to 20°
- Packing unit : 4ea



Locator® Extended Replacement Male

- Retention force : Approx. 6N
- Placement angle compensated up to 20~40°
- Packing unit : 4ea

- Retention force : Approx. 12N
- Placement angle compensated up to 20~40°
- Packing unit : 4ea



Locator® Abutment Components

Locator® Black Processing Male

- Male used only in prosthesis fabrication process
- Packing unit : 4ea



LBPS

Locator® Block Out Spacers

- Used for sealing of the space between the abutment and the denture cap when attaching the overdenture and denture cap in the oral cavity
- Packing unit : 20ea



LBSS

Locator® Impression Coping

- Pick up impression coping for Locator Abutment
- Using closed tray
- Packing unit : Impression coping + Provisional male 1set



LICS

Locator® Lab Analog

- Lab analog for Locator Abutment
- Packing unit : 4ea



LAL50S

Locator® Core Tool

- Used for placing and removing the replacement male in the denture cap
- Separated into three pieces and used as a hand driver for Locator Abutment



LCCT

Locator® Torque Driver

- Torque driver for Locator Abutment





TWLDSK

TWLDLK

O-ring Abutment

- Abutment for overdenture using O-ring attachment
- Placement angle compensated up to 20°
- Tightened with a dedicated outer driver (code : AORD)
- Recommended tightening torque: 30Ncm

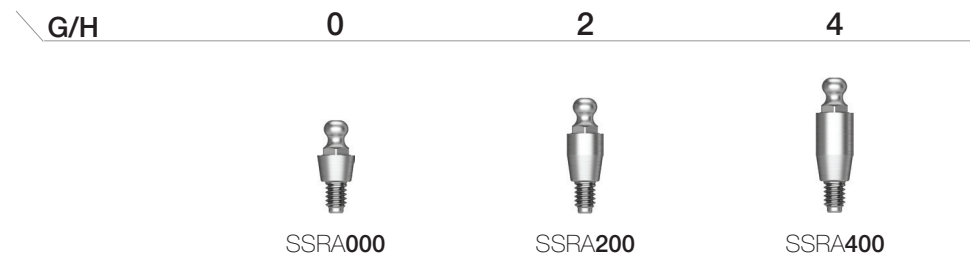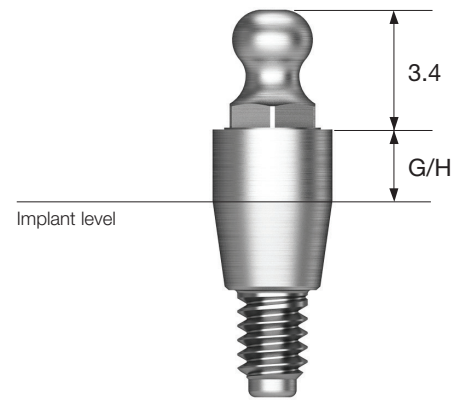

O-ring Abutment Components

O-ring Retainer Cap Set

- O-ring attachment for O-ring Abutment
- O-ring replaced in a metal housing for use
- Packing unit : retainer cap + o-ring



RCS01

O-ring Retainer Set

- Used when vertical dimension is shorter than the retainer cap
- Packing unit : retainer cap + o-ring



RS01

O-ring Set

- O-ring set
- Packing unit : o-ring 5ea



OAON01S

O-ring Lab Analog




- Lab analog for O-ring Abutment



OAL

OneSeal

- Disposable medical devices for internal filling of abutment
- Cut to desired length for use (medical silicone)
- Packing unit : 5ea
- TS Mini : TSSE2250S
- TS Regular : TSSE2350S
- SS Regular : SSSE2650S

H \ D	Ø2.20	Ø2.35	Ø2.65
Type		Long	
			
50	TSSE2250S	TSSE2350S	SSSE2650S



OSSTEM[®]
IMPLANT